



**THE SAFE SOLUTION FOR  
YOUR SMALLEST PATIENTS.**



[www.Sun-Med.com](http://www.Sun-Med.com) | 800-433-2797 | [info@Sun-Med.com](mailto:info@Sun-Med.com)





## DELIVER CONSISTENT PIP AND PEEP WITH EVERY BREATH.

The SafeT T-Piece Resuscitator is the safe solution for emergency neonatal resuscitation and CPAP. Delivering continuous oxygen flow alongside pre-set PIP and PEEP pressures enables clinicians to safely provide the necessary ventilatory support during resuscitation of infants. Designed specifically for use on the delicate lungs of premature babies and newborns under 10 kg/22 lbs.

Safely ventilate  
the fragile lungs  
of neonates.

**SafeT T-Piece Evaluations Available**

Request a sample or product demonstration  
for your NICU + labor/delivery team today.

Provide  
consistent  
ventilatory  
pressure.

Pre-set  
PIP and PEEP  
settings prior  
to use.



Integrated, color-coded pressure  
manometer for effective monitoring  
of ventilatory pressures.




PEEP dial allows clinician to set PEEP;  
ventilation is provided by occluding  
PEEP dial hole intermittently at the  
desired respiratory rate allowing for  
inspiration and exhalation to occur.



Use optional tethered override cap  
on inspiratory pressure controller  
when pressures greater than  
40 cm H<sub>2</sub>O are needed.

# SIDE-BY-SIDE

Current NRP guidelines recommend the use of either self-inflating, flow-inflating or T-piece resuscitators for manual ventilation of neonates. This side-by-side comparison highlights the additional features of a T-piece resuscitator compared to a BVM and hyperinflation system.

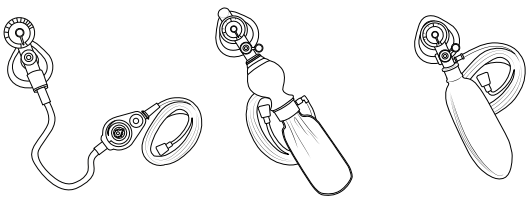

| Neonatal Resuscitator Features + Functionality         | SafeT | BVM | Hyperinflation |
|--|-------|-----|----------------|
| Continuous O <sub>2</sub> Flow                         | ●     |     | ●              |
| Pre-Set + Deliver Consistent PIP                       | ●     |     |                |
| Pre-Set + Deliver Consistent PEEP                      | ●     | ●   | ●              |
| Can Provide CPAP                                       | ●     |     | ●              |
| Pre-set Maximum Pressure + Adjustability During Resus. | ●     |     |                |
| Integrated Pressure Manometer                          | ●     | ●   | ●              |
| Pressure Relief Valve                                  | ●     | ●   | ●              |
| Hand Compression Required                              |       | ●   | ●              |
| Portable   | ●     | ●   | ●              |

## CONFIGURATIONS



## SAFET T-PIECE RESUSCITATOR

| ITEM  | DESCRIPTION                           | MASK SIZE            | PK |
|-------|---------------------------------------|----------------------|----|
| ST700 | SafeT T-Piece Resuscitator            | X-Small Neonate mask | 10 |
| ST701 | SafeT T-Piece Resuscitator            | Small Neonate mask   | 10 |
| ST702 | SafeT T-Piece Resuscitator            | Small Infant mask    | 10 |
| ST703 | SafeT T-Piece Resuscitator            | Infant mask          | 10 |
| ST704 | SafeT T-Piece Resuscitator            | Neonate Round mask   | 10 |
| ST705 | SafeT T-Piece Resuscitator            | Infant Round mask    | 10 |
| ST707 | SafeT T-Piece Resuscitator            | No Mask              | 10 |
| ST710 | SafeT T-Piece Resuscitator w/ Pop Off | X-Small Neonate mask | 10 |
| ST711 | SafeT T-Piece Resuscitator w/ Pop Off | Small Neonate mask   | 10 |
| ST712 | SafeT T-Piece Resuscitator w/ Pop Off | Small Infant mask    | 10 |
| ST713 | SafeT T-Piece Resuscitator w/ Pop Off | Infant mask          | 10 |
| ST714 | SafeT T-Piece Resuscitator w/ Pop Off | Neonate Round mask   | 10 |
| ST715 | SafeT T-Piece Resuscitator w/ Pop Off | Infant Round mask    | 10 |
| ST717 | SafeT T-Piece Resuscitator w/ Pop Off | No Mask              | 10 |

