

safeT™

T-PIECE RESUSCITATOR

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE

The SafeT™ T-Piece Resuscitator is the safe solution for neonatal resuscitation and CPAP. Delivering continuous oxygen flow while providing consistent PIP and PEEP pressures allows the clinicians to provide the necessary ventilatory support during resuscitation. SafeT helps prevent over-pressurization that can lead to barotrauma and pneumothorax in neonates.

Integrated, color-coded pressure manometer for effective monitoring of ventilatory pressures

PEEP dial allows clinician to set PEEP; ventilation is provided by occluding PEEP dial hole intermittently at the desired respiratory rate allowing for inspiration and exhalation to occur



Variety of mask sizes available to meet patient and hospital needs

Standard pressure relief valve designed to lower the risk of over-pressurization


Oxygen tubing with red Fits-All universal connector

safeT™ T-PIECE RESUSCITATOR

- Allows clinician to deliver consistent PIP and PEEP pressures
- Consistent PIP and PEEP ensures greater safety when ventilating babies' fragile lungs, particularly in an emergency situation
- Integrated pressure manometer provides clear visualization of ventilatory pressure ranges
- Pressure relief valve helps prevent over-pressurization of lungs
- Able to provide CPAP
- Includes a cushioned face mask for proper sealing during mask ventilation
- Easy to read controller and manometer makes set up simple and quick
- Universal red Fits-All O₂ connector allows immediate visualization of proper oxygen connection

Expandable circuit that allows up to 18" in length

Inspiratory pressure controller allows clinician to set the PIP

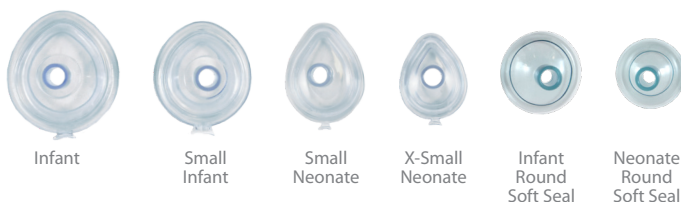

ITEM	DESCRIPTION	MASK	PK
ST700	SafeT T-Piece Resuscitator	X-Small Neonate	10
ST701	SafeT T-Piece Resuscitator	Small Neonate	10
ST702	SafeT T-Piece Resuscitator	Small Infant	10
ST703	SafeT T-Piece Resuscitator	Infant	10
ST704	SafeT T-Piece Resuscitator	Neonate Round	10
ST705	SafeT T-Piece Resuscitator	Infant Round	10
ST707	SafeT T-Piece Resuscitator	—	10